


PharmaPort Gloveports, Gloves and Equipment Modification

Extract Technology can now supply its PharmaPort gloveport as a low cost component part to any Isolator, Glovebox or similar equipment. They are also able to offer a full up-grade package which will allow the Extract team to modify any existing equipment to incorporate the ergonomically designed gloveport into their process.

Providing our latest ergonomic PharmaPort gloveports is one such Isolator upgrade, allowing the operators greater freedom of movement due to their oval shape. They are commonly mounted in a 25mm thick clear acrylic panel with a protective film of Hardcoat™ applied for improved scratch and solvent resistance, but can also be supplied in standard toughened laminated glass.

Features include:

- Helps eradicate contamination hang-up around the operator access glove and the window panel
- Improved operator safety through clamping the glove to the port
- Anti static coating on the acrylic screen
- Improved ergonomic design
- Range of styles and types of gloves to suit all requirements
- Safe change facility for replacing the gloves
- Can be installed into any equipment
- Smooth cleanable internal surface.

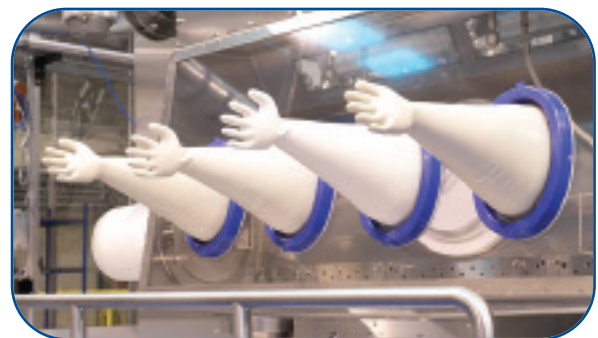

If you are an existing Extract customer, then we will have the original drawing layout for the Isolators we have supplied and can therefore upgrade them very simply and easily.

Even if the equipment has been supplied by another company, a simple site visit to measure the components will generally suffice, allowing us to provide a quotation.

Glove/Gauntlet Supply:

Extract Technology can offer an extensive range of gloves and gauntlets in a variety of materials. These include Natural Latex, Neoprene, Hypalon, Butyl, Polyurethane and Viton covered applications.

All gauntlets/gloves are available with a number of port diameters, thicknesses, lengths and hand sizes to suit individual operator requirements. Extract can now also offer bespoke 3 part gauntlets which give operators sleeve protection whilst allowing the flexibility of a disposable hand glove which allows for precise and intricate operations.

We can also operate a system whereby scheduled deliveries are made to customers, therefore guaranteeing an uninterrupted supply of gloves/gauntlets in situations where continual production is essential.



Glove Integrity Test Kit

To provide glove integrity, a purpose designed Extract Technology Glove Test Kit is also available. This comprises of a machined cast nylon port blanking plug which would include a pneumatic seal to suit the internal bore of the glove port. Mounted to the blanking plate would be a compressed air supply point and Dwyer magnehelic gauge for visual indication of the glove condition.

